




# Instructions for set up of the Starlock cushion




Before you start please ensure the blue valve is in line with the black pipe then press the button for approximately 10 secs

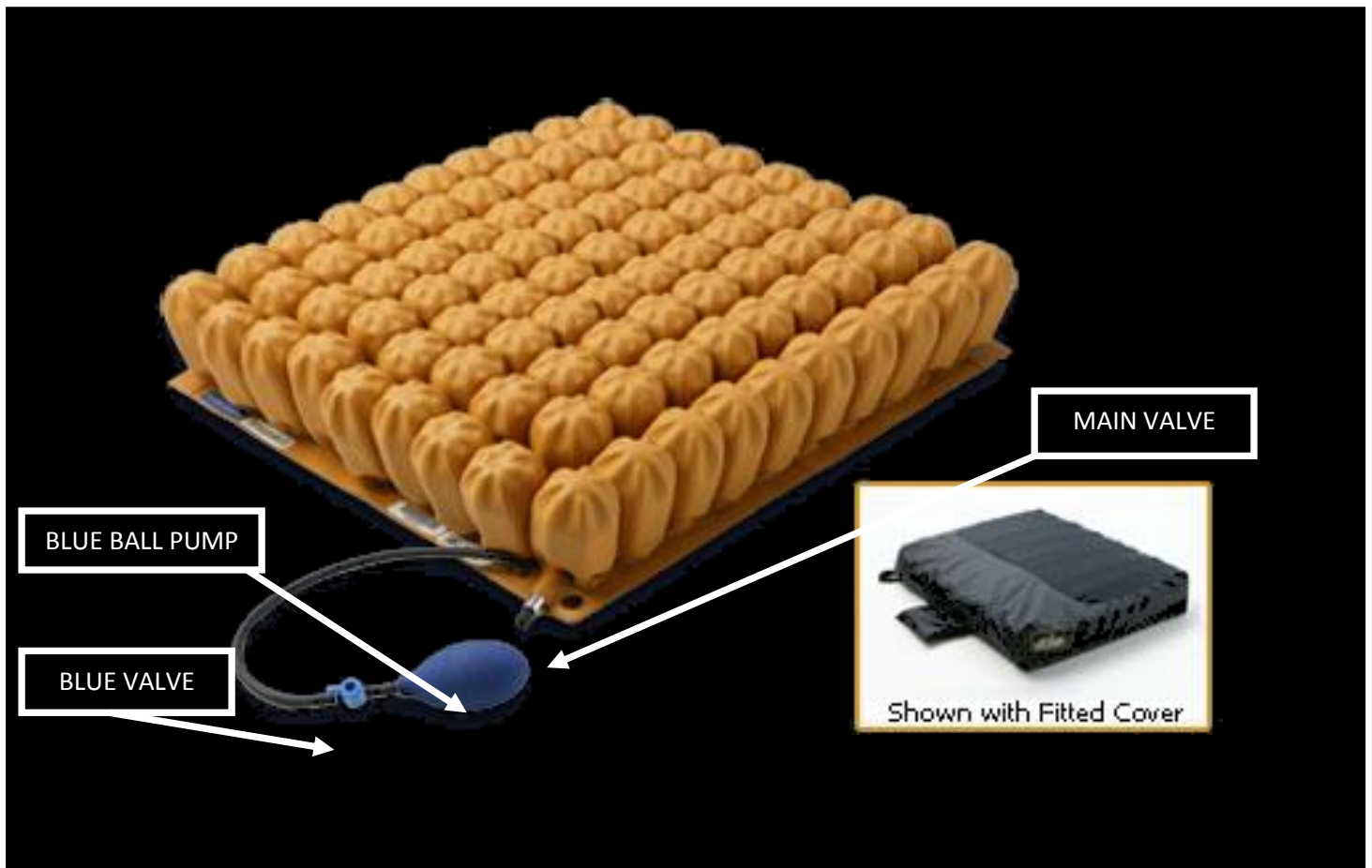


## INITIAL SET-UP

1. Undo the main valve and pump air into the cushion using the pump provided, continue until the cushion is over inflated. Approximately 25 pumps.

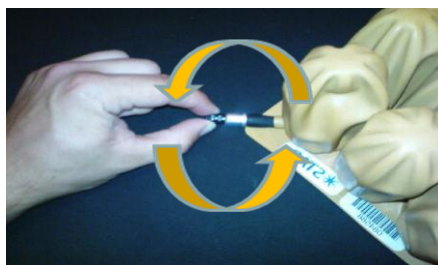

2. Transfer or sit your client onto the over inflated Starlock
3. Once comfortable place your hand under your client specifying their most prominent area whilst seated. (Usually one of the IT's)

**Place hand in here  
between client and  
cushion**

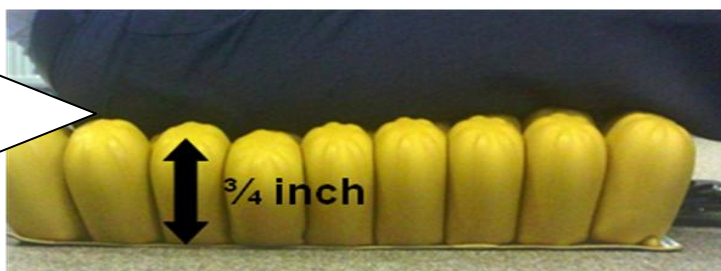

4. Whilst in the position undo the main valve anti-clockwise and release air out of the cushion until there is  $\frac{3}{4}$  inch between your clients prominence and the base of the cushion.

5. Close off valve- clockwise




The Starlock can be personalised to suit the client in 2 key ways:

### **OPTION 1**

If your client has a postural need or just requires stabilisation follow these next 5 steps


#### **1. Position your client in the desired position**

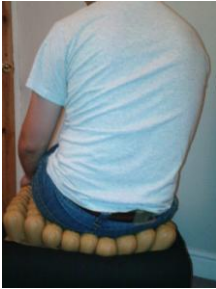

#### **2. Pump the Blue ball pump between 12-15 times**



#### **3. Twist blue valve horizontally**

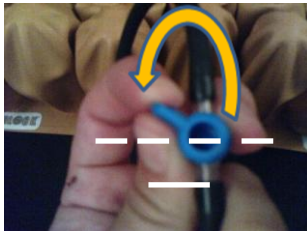

**4. If client ISN'T in a good position repeat steps 1-3 until desired result achieved.**

**5. Once the desired position has been met, secure the Ball pump and pipe into the black pouch at the front of the cushion.**

*YOUR CLIENT NOW HAS THE BENEFIT OF HIGH RISK PRESSURE REDUCTION  
SUPPORTED WITH STABLE POSITIONING*


## OPTION 2

If your client has a Vulnerable Area that needs pressure taking off please follow these next 6 steps:

1. Take pressure off the vulnerable area by rolling up one of the below:  
Pillow Case, Tea Towel, Hanky, Tissue Paper
2. Place this between the clients vulnerable area and the cushion, this may be quite sensitive but it will only take a few seconds
3. Pump the blue ball pump between 12-15 times



4. Turn the blue valve horizontally across the pipe

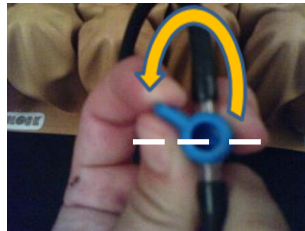

5. Take out the object you used under the vulnerable area

6. Having removed all contact in the key area and deflated pressure, store the ball pump along with the pipe into the pouch



*YOUR CLIENT NOW HAS THE BENEFIT OF HIGH RISK PRESSURE REDUCTION WITH ADDED BENEFIT OF NO PRESSURE ON THAT OPEN WOUND OR VULNERABLE AREA*

If you have any queries please contact your Area Sales Manager  
or Head office on +44 (0)1531 635678